



Macrosep® Advance Centrifugal Devices

- Concentrate and purify samples up to 20 mL.
- Provides recoveries typically > 90%.
- Built-in deadstop prevents spinning to dryness.

Ordering Information

Macrosep Advance Centrifugal Devices

| Description | 6/pkg | 24/pkg | 100/pkg |
|--------------------------|--------------|---------------|----------------|
| 1K Omega™, Yellow | MAP001C36 | MAP001C37 | MAP001C38 |
| 3K Omega, Gray | MAP003C36 | MAP003C37 | MAP003C38 |
| 10K Omega, Blue | MAP010C36 | MAP010C37 | MAP010C38 |
| 30K Omega, Red | MAP030C36 | MAP030C37 | MAP030C38 |
| 100K Omega, Clear | MAP100C36 | MAP100C37 | MAP100C38 |
| 0.2 µm Supor®, Aqua | — | MAPM02C67 | MAPM02C68 |
| 0.45 µm Supor, Wildberry | — | MAPM45C67 | MAPM45C68 |

Table of Contents

| Section | Page |
|--|------|
| Introduction | 1 |
| Filtration Principles for Macrosep Advance Centrifugal Devices | 1 |
| Specifications | 2 |
| Applications | 3 |
| Choosing the Appropriate Macrosep Advance Centrifugal Device | 4 |
| Components | 6 |
| Macrosep Advance Centrifugal Device Operation | 7 |
| Instructions for Use | 7 |
| Pre-Rinsing (Optional) | 7 |
| Non-Specific Adsorption | 8 |
| Diafiltration (Desalting and Buffer Exchange) | 9 |
| Optimization | 9 |
| Complementary Products | 12 |
| Warning | 12 |

Note: *The procedures herein are intended only as a guide. Users should always verify product performance with their specific applications under actual use conditions. If you have questions about the information presented in this guide, please contact our Technical Service department.*

Introduction

Macrosep Advance centrifugal devices provide rapid and efficient concentration and purification of up to 20 mL of biological samples. The unique design maximizes filtration area to process samples quickly while maintaining a gentle concentration environment to preserve protein activity and conformation. The wide selection of ultrafiltration molecular weight cut-off (MWCO) devices incorporate Omega membrane which is very low in protein and nucleic acid binding. Ultrafiltration devices are ideal for concentrating small peptides, oligonucleotides, nucleic acids, enzymes, antibodies, and other similar macromolecules. Macrosep Advance centrifugal filters are also available in 0.2 and 0.45 μm pore sizes containing Pall's Supor polyethersulfone membrane for low protein and nucleic acid binding with high chemical compatibility. The microporous membrane selections are ideal for microorganism concentration, sample clarification, removal of particulates and colloids, and gentle elution of nucleic acid from agarose gels.

Filtration Principles for Macrosep Advance Centrifugal Devices

Centrifugation provides the driving force for filtration. Ultrafiltration devices are typically centrifuged between 1,000 to 5,000 x g. Biomolecules larger than the nominal MWCO of the membrane are retained in the sample reservoir while solutions and low molecular weight molecules pass through the membrane into the filtrate receiver. Microfiltration membrane devices can be centrifuged up to 14,000 x g. Similarly, particulate larger than the membrane pore size are retained in the sample reservoir while solutions and particulate smaller than the pore size pass through into the filtrate receiver.

Specifications

Materials of Construction

Filtration Media:

Ultrafiltration: Omega membrane (modified polyethersulfone)

Microfiltration: Supor membrane (polyethersulfone)

Sample Reservoir, Filtrate Receiver, and Cap: Polypropylene

Paddle: Polyethylene

Effective Filtration Area

7.2 cm² (1.12 in.²)

Dimensions

Diameter: 29 mm (1.2 in.)

Length: 12.0 cm (4.7 in.)

Capacities

Maximum Sample Volume: 20 mL

Maximum Filtrate Receiver Volume: 22 mL

Hold-up Volume: 80 µL

Dead Stop Volume:

34° fixed angle 1.5 mL

45° fixed angle 1.2 mL

Swinging bucket 450 µL

Operating Temperature Range

0-40 °C (32-104 °F)

pH Range

Ultrafiltration: 1-14

Microfiltration: 1-14

Maximum Centrifugal Force

Ultrafiltration: 5,000 x g

Microfiltration: 14,000 x g

Centrifuge

Fits centrifuges that accept standard 50 mL conical-end tubes

Sanitization

Provided non-sterile, may be sanitized by filtering 70% ethanol through the device prior to use.

Applications

Macrosep Advance centrifugal devices with ultrafiltration membrane can be used for:

- Concentrate and desalt proteins and nucleic acids
- Buffer exchange or salt removal of chromatography fractions
- Harvest biomolecules from cell culture media
- Virus concentration or removal
- Crude fractionation of protein mixtures
- Remove debris and particulate from cell lysates

Macrosep Advance centrifugal devices with microfiltration membrane can be used for:

- Separate DNA from agarose gels
- Separate proteins, oligonucleotides, and RNA from polyacrylamide gels
- Clarify samples before HPLC analysis
- Remove cells from media prior to analysis
- Filtration of biological samples
- Collect and wash treated particles or beads
- Fill with a chromatographic medium for analytical procedures or process development

Choosing the Appropriate Macrosep Advance Centrifugal Device for Ultrafiltration Applications

Protein Applications

For maximum retention, select a Macrosep Advance centrifugal device with a MWCO 3 to 6 times less than the molecular weight of the protein to be retained. For example, for a 150K protein, a 30K Macrosep Advance centrifugal device would be the appropriate selection.

DNA Applications

The molecular weight of a strand of DNA can be estimated by multiplying the number of bases by 340 for single stranded DNA, and the number of base pairs by 680 for double stranded DNA. Once the molecular weight of the DNA is estimated, select a Macrosep Advance centrifugal device with a molecular weight cutoff 3 to 6 times less than the molecular weight of the DNA to be retained. For example, to retain a 2 kilobase (Kb) double stranded DNA fragment: $2000 \times 680 = 1,360,000$ Daltons = 1360K Daltons; a 100K Macrosep Advance centrifugal device would be the appropriate selection.

Choosing the Appropriate Macrosep Advance Centrifugal Device for Ultrafiltration Applications *(continued)*

The table below is a guide for initial selection of Macrosep Advance centrifugal devices MWCOs for retention of proteins and nucleic acids. Ionic conditions, molecular conformation, and protein:protein interactions can affect retention of biomolecules. We recommend pretesting retentivity with your biomolecular solution.

Table 1. Macrosep Advance Selection

| Macrosep Advance Device MWCO | Recommended g-force | Biomolecule Molecular Weight or Size | Nucleic Acid Base Pair (ds) | Bases (ss) |
|---|--------------------------------|---|--|-------------------|
| 1K | 3,000-5,000 x g | 3K-10K | 5-16 bp | 9-32 bs |
| 3K | 3,000-5,000 x g | 10K-20K | 16-32 bp | 32-65 bs |
| 10K | 3,000-5,000 x g | 30K-90K | 50-145 bp | 95-285 bs |
| 30K | 3,000-5,000 x g | 90K-180K | 145-285 bp | 285-570 bs |
| 100K | 1,000-3,000 x g | > 300K | 475-1,450 bp | 950-2,900 bs |


Components

Each Macrosep Advance centrifugal device consists of a screw-on cap, sample reservoir containing a paddle with sealed membrane on both sides, and a filtrate receiver tube.

The insert and standard centrifuge tube design provides maximum stability for handling and centrifugation. The filtrate receiver tube provides graduations to measure buffer and samples plus a large area to clearly label sample identification.

Figure 1

Macrosep Advance Centrifugal Device

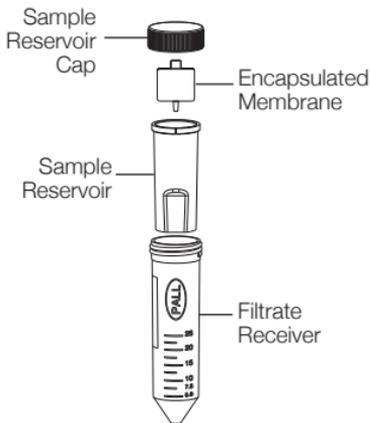

Macrosep Advance Centrifugal Device Operation

Instructions for Use

1. Remove cap and pipette 5 to 20 mL sample into sample reservoir and replace cap to prevent evaporation during centrifugation.
2. Place device into centrifuge that accepts 50 mL conical-end tubes. Always counterbalance the rotor with another Macrosep Advance centrifugal device containing the equivalent sample volume.
3. Spin device at recommended force for required time.
 - a. Ultrafiltration: Spin at 1,000 to 5,000 x g, typically for 30 to 90 minutes, to achieve desired concentrate volume. It is recommended that spin time and g-force be determined for each application.
 - b. Microfiltration: Spin at up to 14,000 x g for 1-3 minutes.
4. Remove the device from the centrifuge and recover target of interest retained in sample reservoir or filtrate receiver tube.
 - a. Target of interest in the sample reservoir: Use pipette to transfer concentrated sample to microcentrifuge tube for storage.
 - b. Target of interest in the filtrate receiver: Remove and discard the sample reservoir and tightly cap the filtrate receiver for storage.

Pre-Rinsing (Optional)

For the majority of applications, Macrosep Advance centrifugal devices can be used without pre-rinsing. However, under certain conditions, it may be preferable to remove trace extractables.

Microfiltration devices: Contact with some organic solvents may cause materials to leach out from the device components. If these leachables represent potential assay interferences, they may be removed by filtering 20 mL of the solvent to be used in the application at 14,000 x g for 1 minute. Discard filtrate and repeat.

Ultrafiltration devices: Omega membrane contains trace amounts of glycerine and sodium azide. If these chemicals interfere with an assay, they may be removed by filtering 20 mL deionized water or buffer through the membrane and repeat. If further flushing is required, start with 0.05N NaOH and repeat this procedure. Use the device within 20 minutes to prevent irreversible membrane damage due to dehydration.

Macrosep Advance Centrifugal Device Operation *(continued)*

Non-Specific Adsorption

Omega membranes are made from polyethersulfone specifically modified to minimize protein binding. These membranes provide equivalent or higher recoveries than comparable regenerated cellulose membranes and offer exceptional biological and chemical resistance.

Adsorption to device components is of particular concern when purifying microgram or nanogram levels of protein. Even with the advanced plastics used in Macrosep Advance centrifugal devices, some adsorption may occur with particularly “sticky” proteins and biomolecules. Pre-treating Macrosep Advance centrifugal devices may further reduce non-specific adsorption to the device.

1. Fill sample reservoir with 20 mL of 10% glycerine.
2. Soak overnight at room temperature.
3. Rinse the device with deionized water.
4. Fill the sample reservoir with 20 mL of deionized water and spin. Repeat.

Macrosep Advance Centrifugal Device Operation *(continued)*

Diafiltration (Desalting and Buffer Exchange)

For salt removal or buffer exchange:

1. Concentrate the sample at least tenfold (e.g., 20 mL concentrated to 2 mL).
2. Reconstitute with exchange buffer and reconcentrate tenfold.
3. Repeat this procedure 3 to 5 times to remove 95 to 99% of salt or buffer.

Optimization

Factors Affecting Performance

Variations in flow rates and recovery can be caused by the following: protein concentration (Macrosep Advance centrifugal devices perform optimally at 1 mg/mL or less protein); temperature (slower flow rates occur at colder temperatures); protein:protein interactions that may cause retention of molecules that would normally pass through the membrane; ionic conditions; and size or conformation of the molecule.

Macrosep Advance Centrifugal Device Operation *(continued)*

Table 2. Effects of Centrifugal Force on Concentration Times

| MWCO | Solute | Time to 25x Concentration (min) | | |
|-------------|---------------------------|--|------------------|------------------|
| | | 1,000 x g | 3,000 x g | 5,000 x g |
| 1K | Ubiquitin (0.25 mg/mL) | >180 min | >180 min | 165 min |
| 3K | Cytochrome C (0.25 mg/mL) | 720 | 300 | 180 |
| 10K | Albumin (1 mg/mL) | 150 | 60 | 45 |
| 30K | IgG (1 mg/mL) | 90 | 60 | 45 |
| 100K | Thyroglobulin (1 mg/mL) | 90 | 60 | 30 |

Macrosep Advance Centrifugal Device Operation *(continued)***Table 3.** Effects of Starting Protein Concentration

| MWCO | Solute | Time to 25x Concentration (min) |
|-------------|--------------------------|--|
| 1K | Ubiquitin 1 mg/mL | 105 |
| | Ubiquitin 0.5 mg/mL | 105 |
| | Ubiquitin 0.1 mg/mL | 105 |
| 3K | Cytochrome C 1 mg/mL | 210 |
| | Cytochrome C 0.5 mg /mL | 180 |
| | Cytochrome C 0.1 mg/mL | 120 |
| 10K | Albumin 1 mg/mL | 45 |
| | Albumin 0.5 mg /mL | 45 |
| | Albumin 0.1 mg/mL | 45 |
| 30K | IgG 1 mg/mL | 45 |
| | IgG 0.5 mg /mL | 30 |
| | IgG 0.1 mg/mL | 30 |
| 100K | Thyroglobulin 1 mg/mL | 30 |
| | Thyroglobulin 0.5 mg /mL | 30 |
| | Thyroglobulin 0.1 mg/mL | 30 |

Complementary Products

- Pall Life Sciences offers centrifugal devices for processing the following sample volumes:

| Device | Sample Volume |
|--------------------------|---------------|
| Nanosep® Device | up to 0.5 mL |
| Microsep™ Advance Device | up to 5 mL |
| Macrosep Advance Device | up to 20 mL |
| Jumbosep™ Device | up to 60 mL |

- **Minimate™ Tangential Flow Filtration Devices** are typically used for the concentration or diafiltration of 100 mL to 5 liter samples.
- **BioTrace™ and Biodyne® Transfer Membranes** offer precise performance and compatibility with nearly every detection system available.
- **AcroPrep™ and AcroPrep Advance 96-well Filter Plates with Supor and Omega Membranes** exhibit low binding capacities for protein and nucleic acid purification.
- **Filtration Devices with Supor Membrane** are sterile, ready-to-use, and maximize sample recoveries with low protein-binding membrane and low hold-up volumes.

WARNING

Employment of the products in applications not specified, or failure to follow all instructions contained in this product information insert, may result in improper functioning of the product, personal injury, or damage to property or the product. See Statement of Warranty in our most recent catalog.

ATTENTION

L'utilisation de nos produits dans des applications pour lesquelles ils ne sont pas spécifiés ou le non-respect du mode d'emploi qui figure sur ce document, peut entraîner un dysfonctionnement du produit, endommager le produit ou d'autres biens matériels ou représenter un risque pour l'utilisateur. Se référer à la clause de garantie de notre catalogue le plus récent.

ACHTUNG

Der Einsatz dieses Produktes in Anwendungen für die es nicht spezifiziert ist, oder das Nichtbeachten einiger, in dieser Bedienungsanleitung gegebenen Hinweise kann zu einem schlechteren Ergebnis, oder Zerstörung des Produktes oder anderer Dinge oder gar zu Verletzungen führen. Beachten Sie auch unsere Garantiebedingungen im aktuellen Katalog.

ADVERTENCIA

El uso de este producto en aplicaciones no especificadas o el no considerar las instrucciones indicadas en la hoja de información del producto puede ocasionar un mal funcionamiento del producto, daños en las instalaciones o en el producto y riesgo para el personal del laboratorio. Consulte el apartado de Garantía en nuestro último catálogo.

ATTENZIONE

L'impiego dei prodotti in applicazioni non specificate, o il mancato rispetto di tutte le istruzioni contenute nel presente bollettino tecnico, potrebbero portare ad un utilizzo improprio del prodotto, ferire gli operatori, o danneggiare le caratteristiche del prodotto stesso. Consultare la dichiarazione di garanzia pubblicata nel nostro più recente catalogo.

警告

当製品情報に記載されていないアプリケーションにおいて製品を使用した場合、あるいは当製品情報に記載されている使用方法に従わない場合は、製品の機能上の不具合、人体への危害、あるいはお客様の財産や製品への損害をまねく恐れがあります。必ず、最新のカタログに記載してある保証約款をご覧ください。



Pall Life Sciences

25 Harbor Park Drive
Port Washington, NY 11050 USA

1.800.521.1520 USA and Canada
(+)800.PALL.LIFE Outside USA and Canada

Visit us on the Web at www.pall.com/lab
E-mail us at LabCustomerSupport@pall.com



**ENABLING A
GREENER
FUTURESM**

To see how Pall is helping enable a greener, safer and more sustainable future, visit www.pall.com/green.

International Offices

Pall Corporation has offices and plants throughout the world in; Argentina, Australia, Austria, Belgium, Brazil, Canada, China, France, Germany, India, Indonesia, Ireland, Italy, Japan, Korea, Malaysia, New Zealand, Norway, Philippines, Poland, Russia, Singapore, South Africa, Spain, Sweden, Switzerland, Taiwan, Thailand, United Kingdom, and Vietnam. Distributors in all major industrial areas of the world. To locate the Pall office or distributor nearest you, visit www.pall.com/contact.

The information provided in this literature was reviewed for accuracy at the time of publication. Product data may be subject to change without notice. For current information consult your local Pall distributor or contact Pall directly.

© 2014, Pall Corporation. Pall, , AcroPrep, Biodyne, BioTrace, Jumbosep, Macrosep, Microsep, Minimate, Nanosep, Omega, and Supor are trademarks of Pall Corporation. The ® indicates a trademark registered in the USA.