

Escherichia coli are common in human bacterial flora. Most of them are harmless, but some known as Shiga Toxic *E. coli* or STEC, can cause severe, sometimes lethal, pathologies.

Main source of infection is contaminated food making it a major concern for food processors.

Among these STEC, several serogroups are more commonly related to human diseases. In Europe, STEC Top 5 serogroups include *E. coli* O157, O26, O103, O111 and O145. In the USA, serogroups O45 and O121 are added to this list making the STEC Top 7.

GeneDisc System Benefits

Rapid — Accelerate the batch release of your short shelf life products and raw materials. While other methods such as immunoassays or culture methods require up to 3 days to results, the Pall GeneDisc method allows a detection of pathogenic STEC in as fast as 10 hours.

Easy to use — GeneDisc solutions are designed for routine use. Implementing PCR (Polymerase Chain Reaction) has never been this easy.

High throughput capability — Process up to 96 samples DNA extractions simultaneously in less than one hour.


Modular — System modularity fits your throughput needs: up to 96 samples can be analyzed in a one hour PCR run.

Choose Your GeneDisc Solutions

Pall GeneDisc Technologies provide you with flexible solutions matching your needs for a reliable monitoring of STEC.

GeneDisc® Technologies

A Flexible Approach for the Detection of Top 5 and Top 7 Shiga Toxic *E. coli*

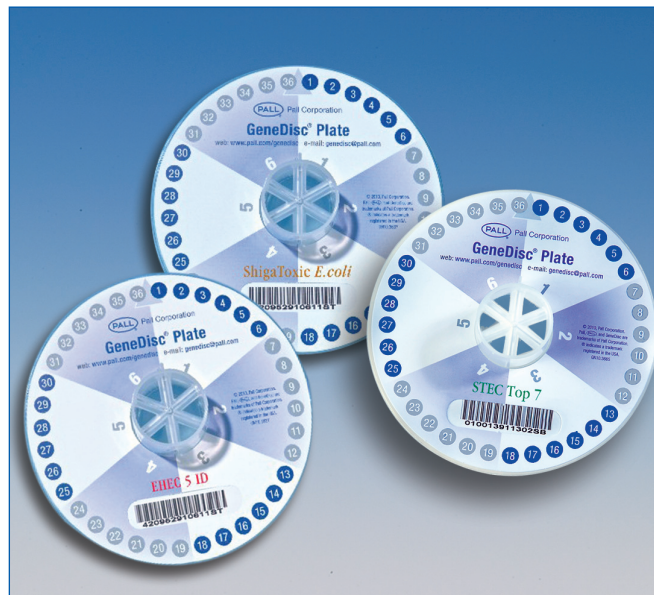

Shiga Toxic *E. coli* Identification

Bacteria	Gram –, motile enterobacteria producing shigatoxins (<i>stx</i>) and intimin (<i>eae</i>)
Food Vehicle	Large variety of food especially beef cattle and dairy products
Disease	Gastroenteritis, Hemorrhagic diarrhea, Hemolytic-uremic syndrome (HUS)
Incidence of Cases (per 100,000 population)	1.90 (FoodNet, 2010)
Notification Rate (per 100,000 population)	0.83 (EFSA, ECDC, 2010)
Related Outbreaks	2 (EFSA, ECDC, 2010) 25 (CDC, 2010)
Related Recalls	8 (FDA & FSIS, 2012)
Related Alerts / Information	22 (RASFF, 2011)

Follow reference method

In line with MLG 5B and ISO/TS 13136 – With this method, a systematic screening based on virulence factors allows to discriminate pathogenic strains from non pathogenic ones. If result is positive, an identification of the Top 7 or Top 5 serogroups is performed.

Test *Salmonella* spp. simultaneously – Analysis with *Salmonella* spp. is available and does not require any additional hands-on time nor enrichment.

Reduce your rate of presumptive positive

Enhanced workflow – This method enables to reduce the number of presumptive positive sample using a cutting edge approach. With this method, all targets – serogroups and virulence factors – are analyzed within one GeneDisc plate.

High level of discrimination – An accurate virulence factor screening based on the association of these factors to serogroups provides a lower rate of presumptive positive than any other available method.



GEN 25/06 – 11/08
ALTERNATIVE ANALYTICAL METHODS FOR AGRIBUSINESS
Certified by AFNOR certification
www.afnor-validation.org



How the System Works



Validations

Method	Part Number	Validation	Matrix	Time to Result (Down to)
STEC	GSTEHEC1XX006	AOAC-PTM	Raw Ground Beef/ Raw Beef Trim	10 h
STEC & <i>Salmonella</i> spp.	GSTECSL2XX006			
STEC Top 7	GTOP7EC106006 GSTECPL206006*			
<i>E. coli</i> O157:H7	GSTEHEC1XX006	NF	Raw Beef Meat	10 h
	GEHECID106006	VALIDATION	Dairy products/ Vegetables products	18 h
	GTOP7EC106006			
<i>E. coli</i> O157:H7 & <i>Salmonella</i> spp.	GSTECSL2XX006		Raw Beef Meat	10 h
	GEHECID106006		Dairy products/ Vegetables products	18 h

All methods are approved using both Extraction Pack Food 1 and Food 2.
*Confirmation tool for the AOAC-PTM certified GeneDisc STEC Top 7 method

Technical Information

Enrichment Time	Down to 8 hours
Sample Preparation Time	< 1 hour for 96 samples
PCR Cycle Time	< 1 hour
Total Turnaround Time	Down to 10 hours
Limit Of Detection	1 bacteria in 25 g of food sample 1 bacteria in 375 g of raw beef meat
Specificity	Wide range of strains tested for inclusivity and exclusivity
Internal Positive Control Per Sample Analysis	Detects presence of inhibitors in each sample DNA extract

Ordering Information

Part Number	Description	Samples/pack
Equipment		
EGDCV3A	GeneDisc Cycler Base Unit	-
EGDSV3A	GeneDisc Cycler Sub Unit	-
EGDUL1A230 (EU)	GeneDisc Ultra-Lyser	-
EGDUL1A120 (US)	OR	
EGDBH96230 (EU)	GeneDisc DryBlock Heater 96	-
EGDBH96120 (US)		
SPSKIT96	Extraction Pack Food 2 Starter Kit	-
Consumables		
PFOOD1100	Extraction Pack Food 1	100
PFOOD2096	Extraction Pack Food 2 (High throughput)	96
GSTEHEC106006	GeneDisc ShigaToxic <i>E. coli</i> *	36
GSTEHEC112006		72
GSTECSL206006	GeneDisc ShigaToxic <i>E. coli</i> *	36
GSTECSL212006	& <i>Salmonella</i> spp.	72
GEHECID106006	GeneDisc EHEC 5 ID (H7, O26, O103, O111, O145)	36
GTOP7EC106006	GeneDisc STEC Top 7 (virulence genes and O157, O26, O103, O111, O145, O45, O121 serogroups)	36
GSTECPL206006	GeneDisc STEC Plus (virulence genes)	36

*Includes identification of pathogenic *E. coli* O157

We also offer a full product range for pathogen detection in food and water and for spoilage organisms in beverage.

Quantitative tests for pathogens in water (*Legionella*, *E. coli*, *Enterococcus*...) are also available.

For more information including part numbers please contact us.



Pall Corporation

Pall GeneDisc Technologies

25 Harbor Park Drive
Port Washington, NY 11050
+1 516 484 3600 telephone
+1 866 905 7255 toll free US

genedisc@pall.com

Visit us on the Web at www.pall.com/genedisc

International Offices

Pall Corporation has offices and plants throughout the world in locations such as: Argentina, Australia, Austria, Belgium, Brazil, Canada, China, France, Germany, India, Indonesia, Ireland, Italy, Japan, Korea, Malaysia, Mexico, the Netherlands, New Zealand, Norway, Poland, Puerto Rico, Russia, Singapore, South Africa, Spain, Sweden, Switzerland, Taiwan, Thailand, the United Kingdom, the United States, and Venezuela. Distributors in all major industrial areas of the world. To locate the Pall office or distributor nearest you, visit www.pall.com/contact.

The information provided in this literature was reviewed for accuracy at the time of publication. Product data may be subject to change without notice. For current information consult your local Pall distributor or contact Pall directly.

© 2014, Pall Corporation. Pall, (PALL), and GeneDisc are trademarks of Pall Corporation. ® indicates a trademark registered in the USA and TM indicates a common law trademark. *Filtration.Separation.Solution.sm* is a service mark of Pall Corporation.